




Patient Teaching

HOME PUMP E(LIPSE

Patient Guidelines:

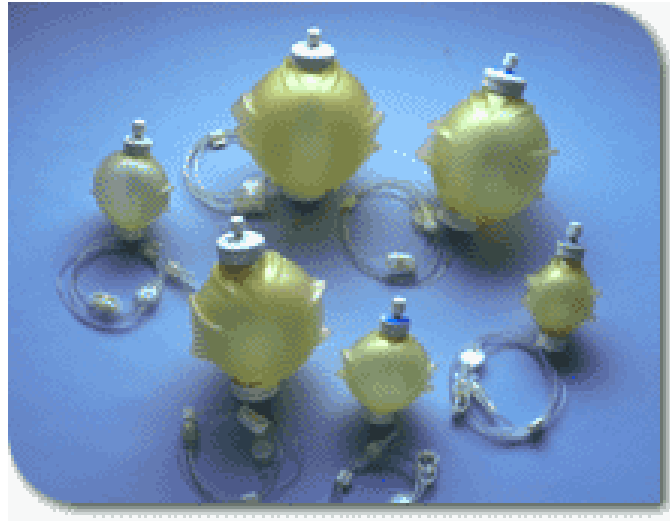


This guide will help you learn the easy steps in using the Eclipse Homepump to administer your medication.

● Storage:

- Your label and your nurse will tell you if your medicine needs to be stored in a refrigerator, Freezer, or at room temperature.
- Your Nurse will tell you when to remove the Homepump from the refrigerator or freezer before use.
 - Remove the Homepump from the refrigerator or freezer _____ hours before use. The colder the fluid the slower it will run.

HOME PUMP ECLIPSE


- Do not use a microwave, oven, or bathwater to warm the fluid. Allow the Homepump to warm at room temperature or hold the Homepump in your hands till desired temperature.

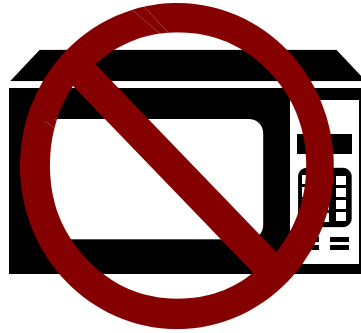

1. Wash your hands.



2. Carefully check the label on the solution bag or bottle for the following before using:
 - a. Your name.
 - b. Name of medication.
 - c. Expiration date.
 - d. Directions for use.



215 Shore Rd
Somers Point NJ
888-646-6379
Dea # 99999999
For Parenteral use only

John Q. Smith Dr Moe, Larry, Cheese
Rx # 99999 Date 12/25/05

Ampicillin 525 mg = 100 ml NSS Eclipse Qty: 11

Directions:
Infuse Ampicillin 525 mg in 100ml NSS over 60 minutes
Every 8 hours via Eclipse Homepump.
Flush before and after each dose using SASH method
Storage : Keep Refrigerated remove from refrigerator
4 hours before use.

Discard after: 1/25/06 RPh JLC
Federal law prohibits dispensing without a
prescription

HOME PUMP E(LIPSE)

3. Prepare a clean workspace by wiping it with alcohol or other anti-germ cleaner.



4. Wash your hands again.




5. Gather your supplies:

2 Alcohol Wipes

2 Saline syringes

1 heparin syringe

Home pump

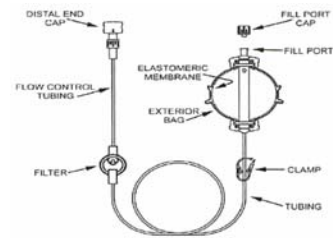


6. Wash your hands again.




HOME PUMP E(LIPSE

7. Wipe the injection cap with alcohol swab on your IV line.

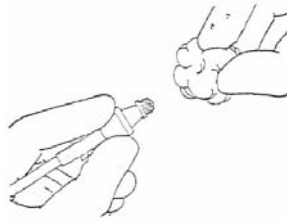


8. Attach a saline syringe to the injection cap and slowly inject saline as you have been instructed.

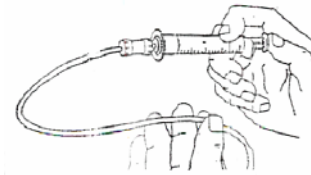


9. Detach syringe from the injection cap.

10. Remove the cover from the end of the IV tubing and connect the tubing to injection cap on the end of your catheter.

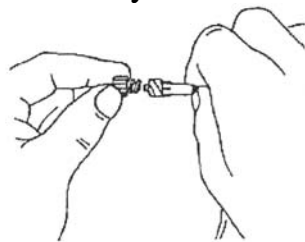


1

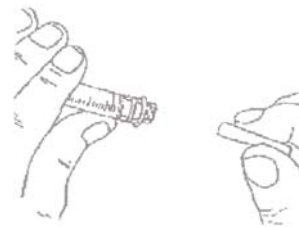

2


11. Open any closed clamps on the Homepump tubing and on your line.

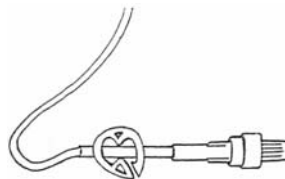



HOME PUMP E(LIPSE

12. Your medicine will start to infuse when you open the clamps.


13. Your infusion is completed when the infusion time has passed and the balloon deflated.

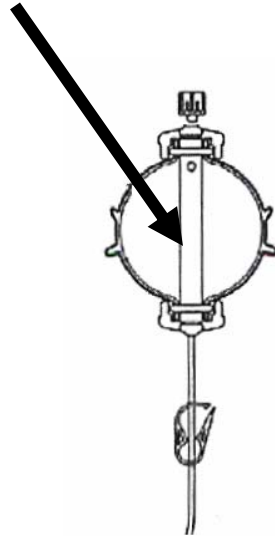

14. Close the clamp on the Home pump and your line remove the end of the tubing from your IV line.

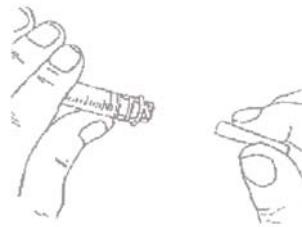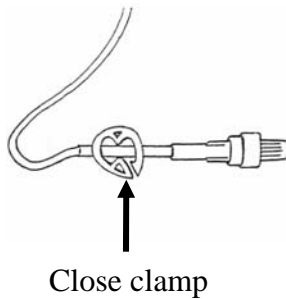


15. Wipe the injection cap with alcohol swab on your IV line.

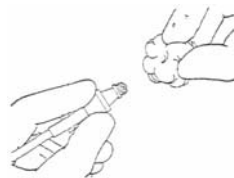


HOME PUMP E(LIPSE

16. Attach a saline syringe to the injection cap and slowly inject saline as you have been instructed. Detach syringe from the injection cap. Clamp your line.

Flush it with _____ ml Saline.

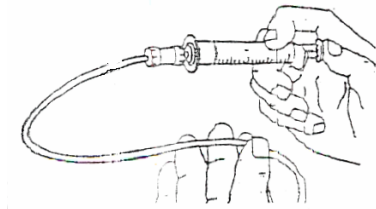

17. Then repeat the same for heparin if indicated. Attach a heparin syringe to the injection cap and slowly inject heparin as you have been instructed. Detach syringe from the injection cap. Flush your line with _____ unit/ml _____ ml of heparin.

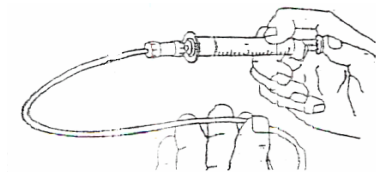

18. Throw away the Homepump ball, and syringes in the regular trash.



Congratulations!

Helpful hints

Remember the word SASH:

S = Saline flush

A = Administer (give) over the number of hours ordered

S = Saline flush after the medication has infused

H = Heparin flush

Always remember if nothing is running through your line clamp it !!!!!

