


Using Dosi-Fuser[®] medicine system

Dosi-fuser[®] is a lightweight, disposable device containing your chemotherapy. We use it to give a slow, continuous infusion into your bloodstream via a central venous catheter or PICC (peripherally inserted central venous catheter). This allows chemotherapy treatment to be given over _____ hours using a Dosi-fuser[®] for home treatment. The chemotherapy is contained in a balloon reservoir, inside a casing. No batteries are needed.



You will also see marks on the bottle that show you when the medicine is full and when it is empty. The tubing is already connected and filled (primed) with the medicine. You will see a small round filter built into the tubing as well as a slide clamp. Towards the end of the tubing is a flow restrictor that needs to be taped to your skin near your IV access.

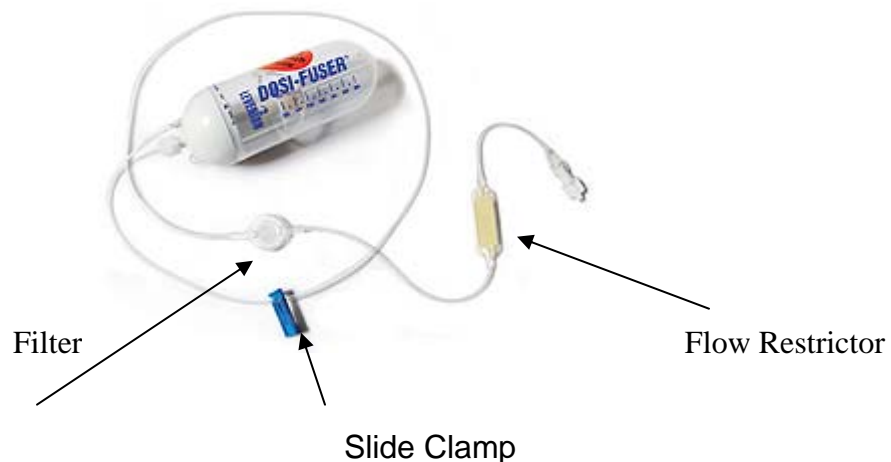

The Dosi-fuser[®] can occasionally stop working or will empty sooner than anticipated. These problems do not happen often but, if they do occur, you should contact us at 1-888-646-6379

You can choose the best way to wear your Dosi-fuser[®]. The Dosi-fuser[®] is supplied with a pouch which you can attach to the inside of a shirt, or jacket, or you may prefer to use the belt/sling.



How Do You Dispose of Your Empty Dosi-fuser[®] ?

Seal the empty infusor with the winged cap connected to prevent leakage of chemotherapy in the yellow bags provided; Home Solutions will pick up your waste with your next delivery. Do not put the infusor in your trash at home.



How Do You Monitor Your Infusor?

You should look at the balloon or scale on the side of the Dosi-fuser[®] every morning to make sure the Dosi-fuser[®] is working correctly.

What To Do If The Dosi-fuser[®] is Faulty or Leaky?

In the unlikely event that your Dosi-fuser[®] is leaking, do not worry; give Home Solutions a call at 1-888-646-6379. You may be asked to follow the steps to detach your Dosi-fuser[®]. Put on gloves and wash any chemotherapy that has spilled onto you or any surface with warm soapy water. If there is some chemotherapy left in the leaking pump, put on a winged cap and seal the Dosi-fuser[®] in a plastic bag, then place it in the yellow plastic bag.

Stopping the infusion

Wash your hands again.



Clean your work area surface with an alcohol prep pad. Then, gather these supplies and place them on your work space:



Gather your supplies:

___ Alcohol Wipes


___ Heparin syringe

___ Saline syringes


Gloves



Close the clamp on your Dosi-fuser[®]

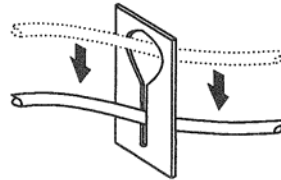

Clamp your line.

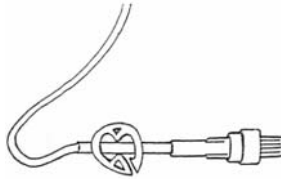


Put on your gloves.



Disconnect the IV tubing from your IV access device


Recap the end.

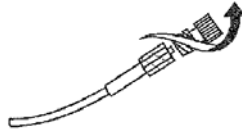


Dispose of the Dosi-fuser[®] in your yellow chemo waste bag and seal.



Wipe the injection cap on your IV line with an alcohol swab for 30 seconds.

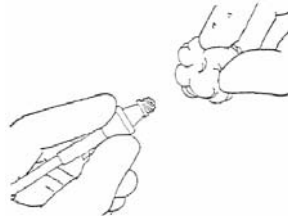

Attach a saline syringe to the injection cap and slowly inject saline as you have been instructed.

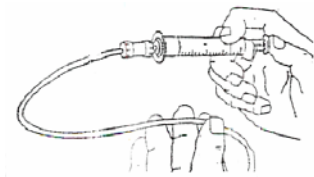

Detach syringe from the injection cap.



Then repeat the same for heparin if indicated. Attach a heparin syringe to the injection cap and slowly inject heparin as you have been instructed. Detach syringe from the injection cap. Flush your line with _____unit/ml _____ml of heparin.

Throw away the Trash, and syringes in the regular trash.



Congratulations!